
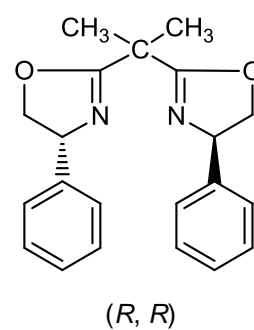
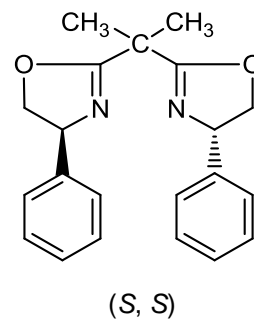


2,2'-イソプロピリデンビス(4-フェニル-2-オキサゾリン)

2,2'-Isopropylidenebis(4-phenyl-2-oxazoline)

E170801P

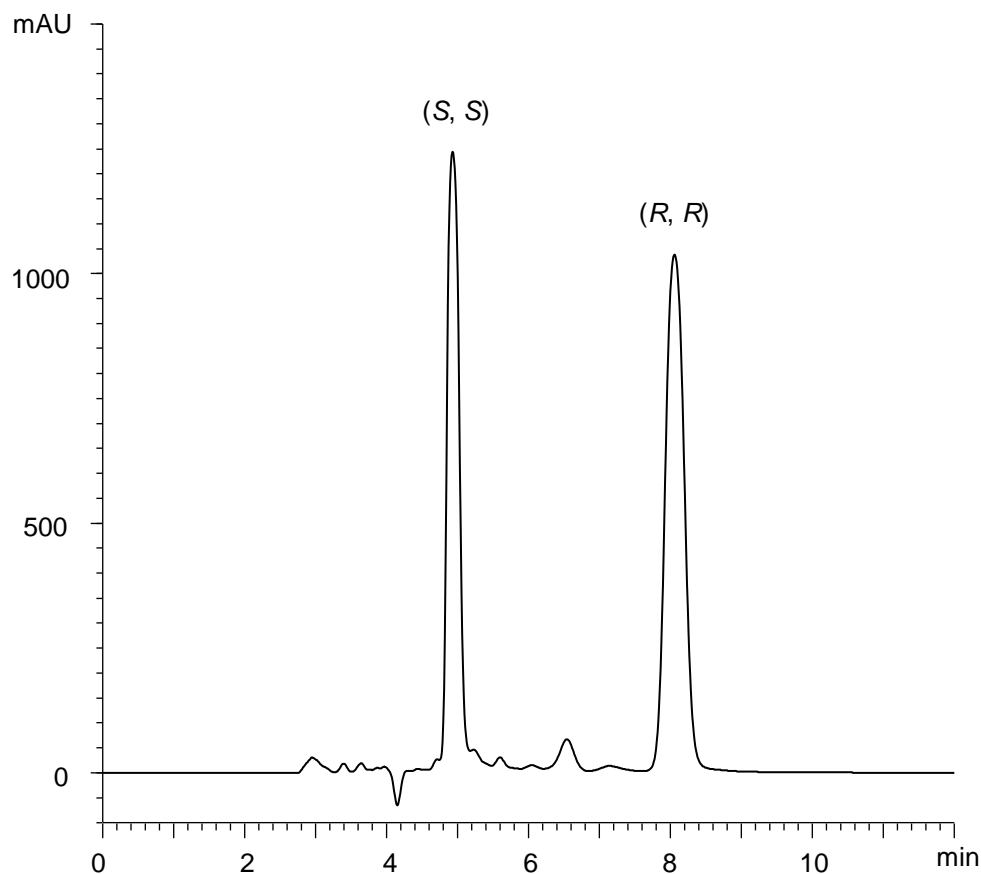

$\alpha$	2.5
$R_s$	7.6

Column : CHIRAL ART Amylose-C (5  $\mu$ m)  
250 X 4.6 mm I.D.

Eluent : *n*-hexane/2-propanol (70/30)

Flow rate : 1.0 mL/min

Temperature : 25°C

Detection : UV at 210 nm

Injection : 10  $\mu$ L (0.5 mg/mL)