


セチリジン (レボセチリジン)


Cetirizine (Levocetirizine)

F180214C

α	1.1
R_s	2.4

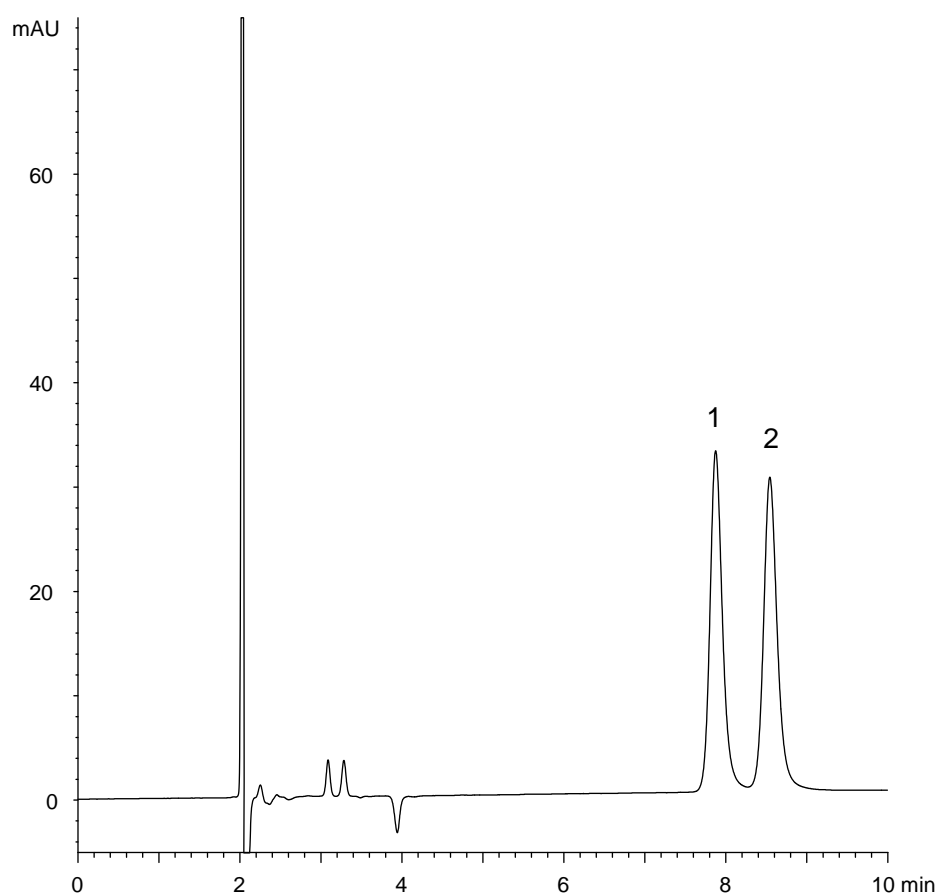


1

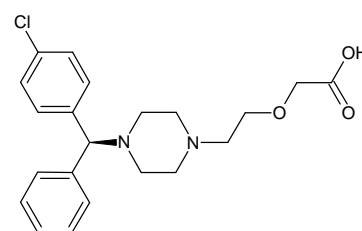

(S)-Cetirizine

2

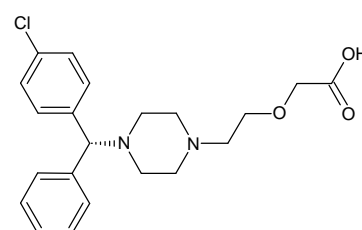

(R)-Cetirizine

(Levocetirizine)

Column : CHIRAL ART Cellulose-SB (3 μ m)
150 X 4.6 mmI.D.
Eluent : acetonitrile/formic acid/diethylamine (100/0.1/0.1)
Flow rate : 1.0 mL/min
Temperature : 25°C
Detection : UV at 230 nm
Injection : 2 μ L (200 μ g/mL)