


モノクローナル抗体
Monoclonal antibody (MAb)

M180604AA

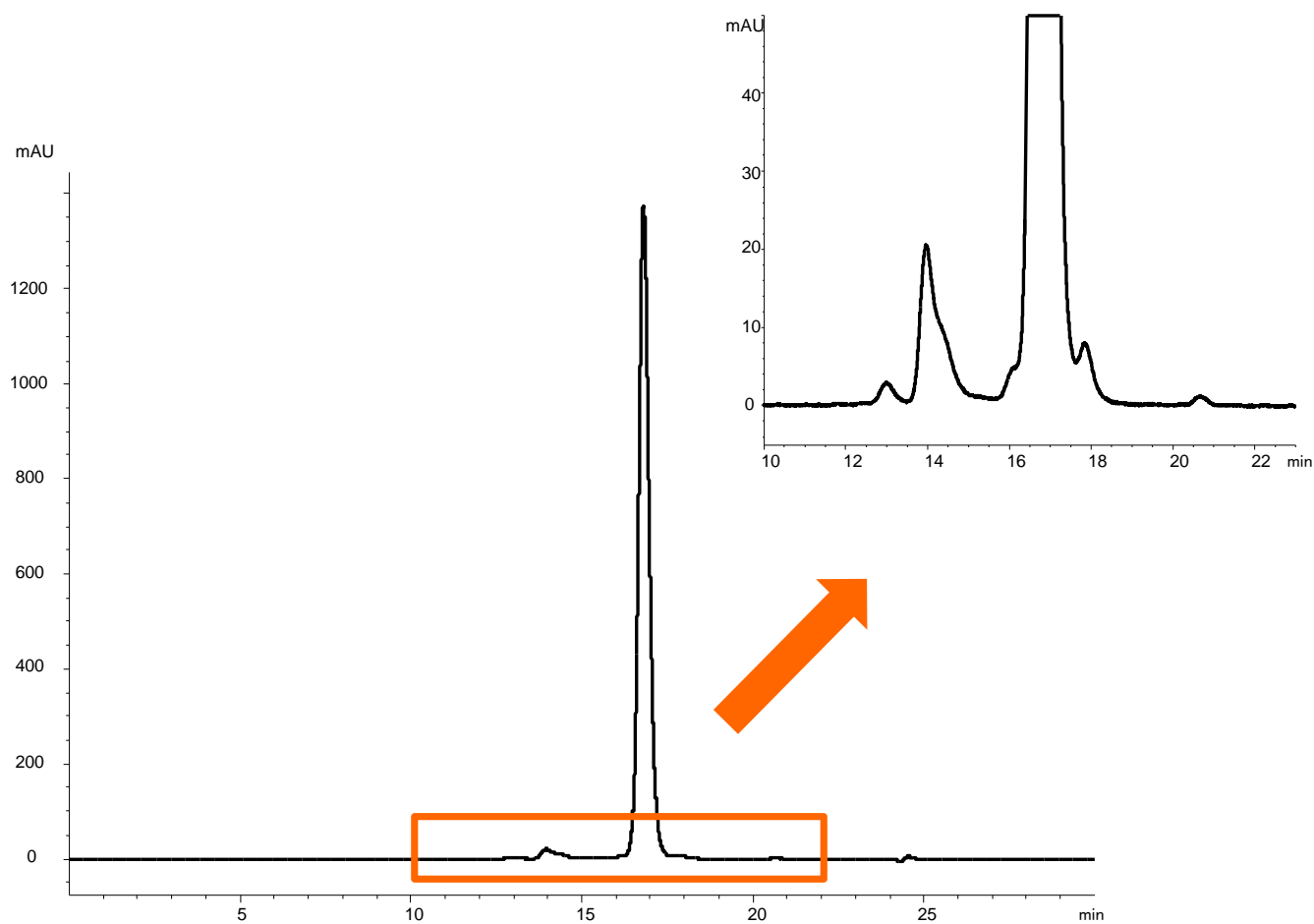

| | |
|-------------|--|
| Column | : YMC-SEC MAB (3 μ m, 25 nm) 300 X 4.6 mmI.D. |
| Eluent | : 0.1 M KH_2PO_4 - K_2HPO_4 (pH 7.0) containing 0.2 M NaCl |
| Flow rate | : 0.165 mL/min |
| Temperature | : 25°C |
| Detection | : UV at 280 nm |
| Injection | : 10 μ L (5 mg/mL Humanized monoclonal IgG1) |