



フェノキシベンズアミン塩酸塩
Phenoxybenzamine Hydrochloride

T180319AP

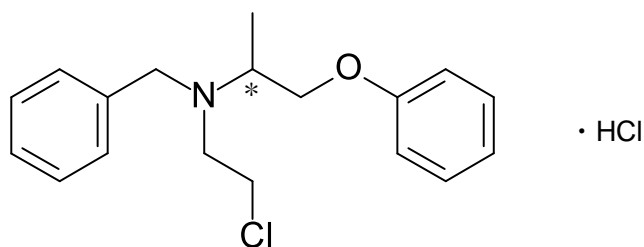

α	1.4
Rs	6.9

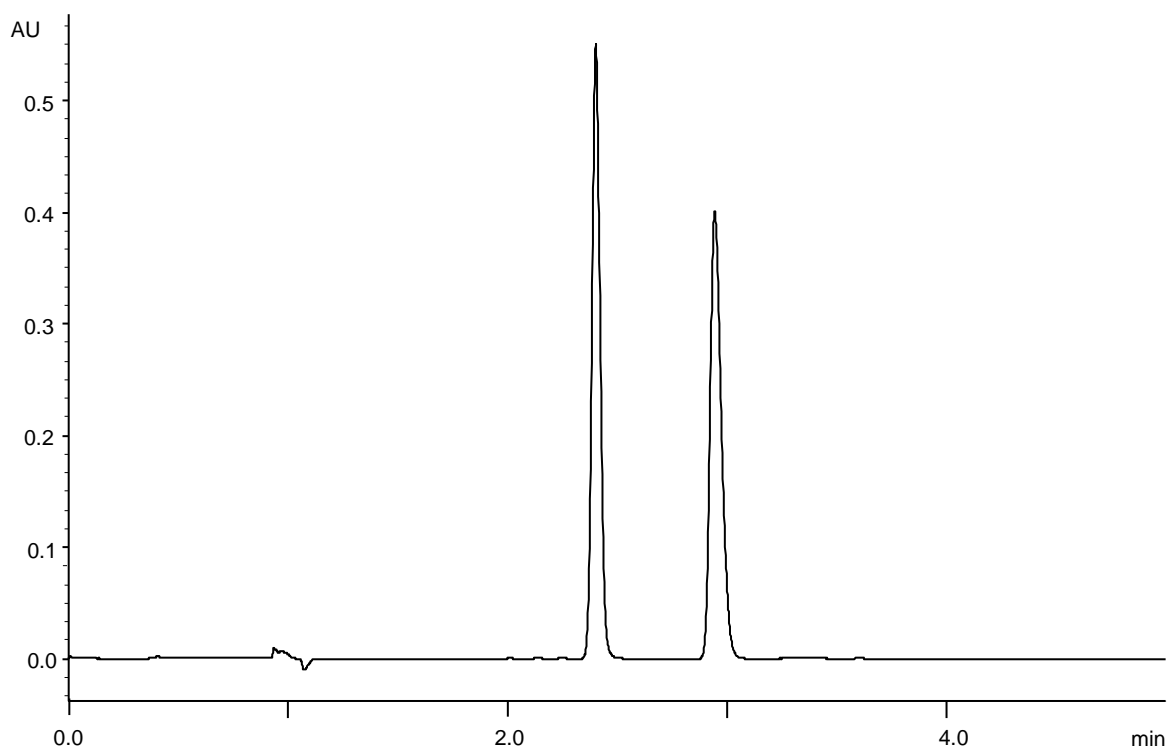

Column : Alcyon SFC CSP Cellulose-SJ (5 μ m)
250 X 4.6 mm I.D.

Eluent : CO₂/2-propanol (85/15)

Flow rate : 3.0 mL/min

Temperature : 35°C

Detection : UV at 210 nm

Backpressure : 13.8 MPa (2000 psi)

Injection : 2 μ L (2 mg/mL)