



インダパミド  
Indapamide

T180319AZ

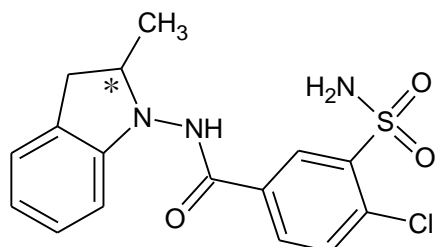

$\alpha$	1.2
$R_s$	3.1

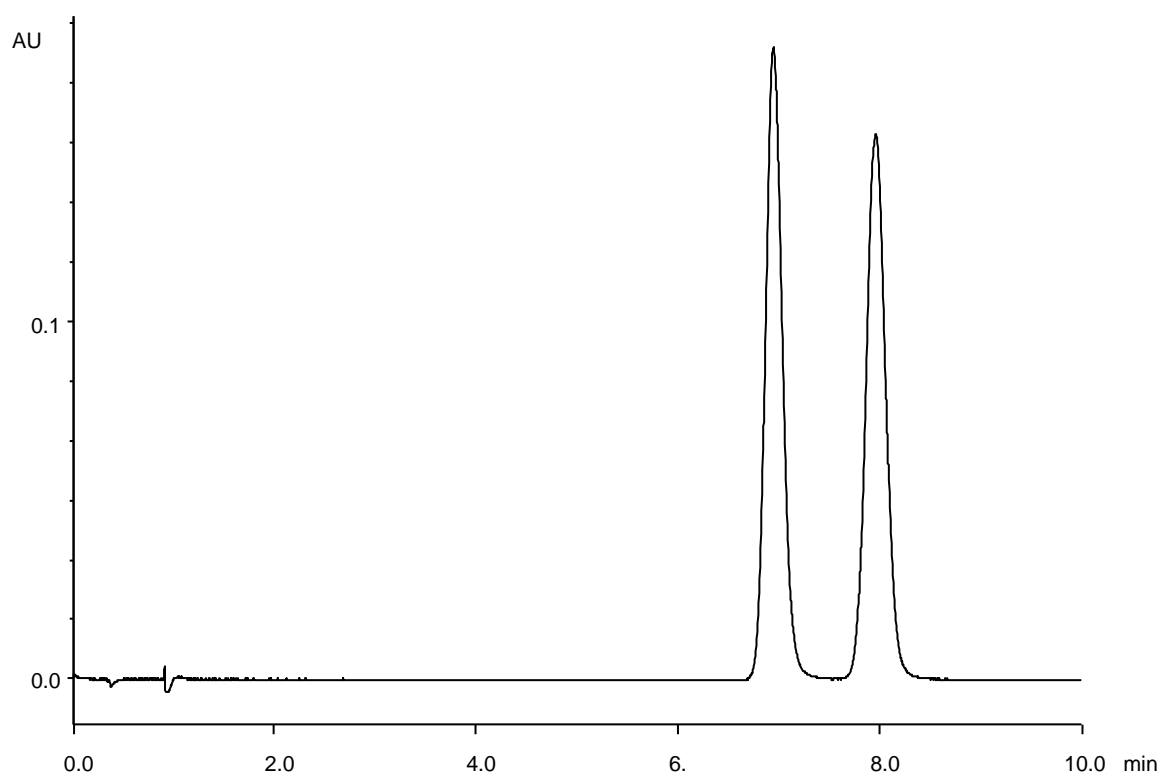

Column : Alcyon SFC CSP Cellulose-SJ (5  $\mu$ m)  
250 X 4.6 mmI.D.

Eluent : CO<sub>2</sub>/ethanol (75/25)

Flow rate : 3.0 mL/min

Temperature : 35°C

Detection : UV at 240 nm

Backpressure : 13.8 MPa (2000 psi)

Injection : 2  $\mu$ L (2 mg/mL)