



1,2-ビス[(2-メトキシフェニル)フェニルホスフィノ]エタン
1,2-Bis[(2-methoxyphenyl)phenylphosphino]ethane


T180515B

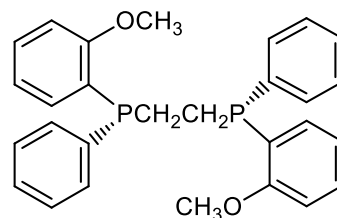

(R, R)

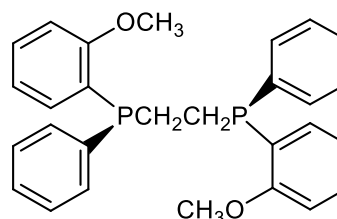

(S, S)

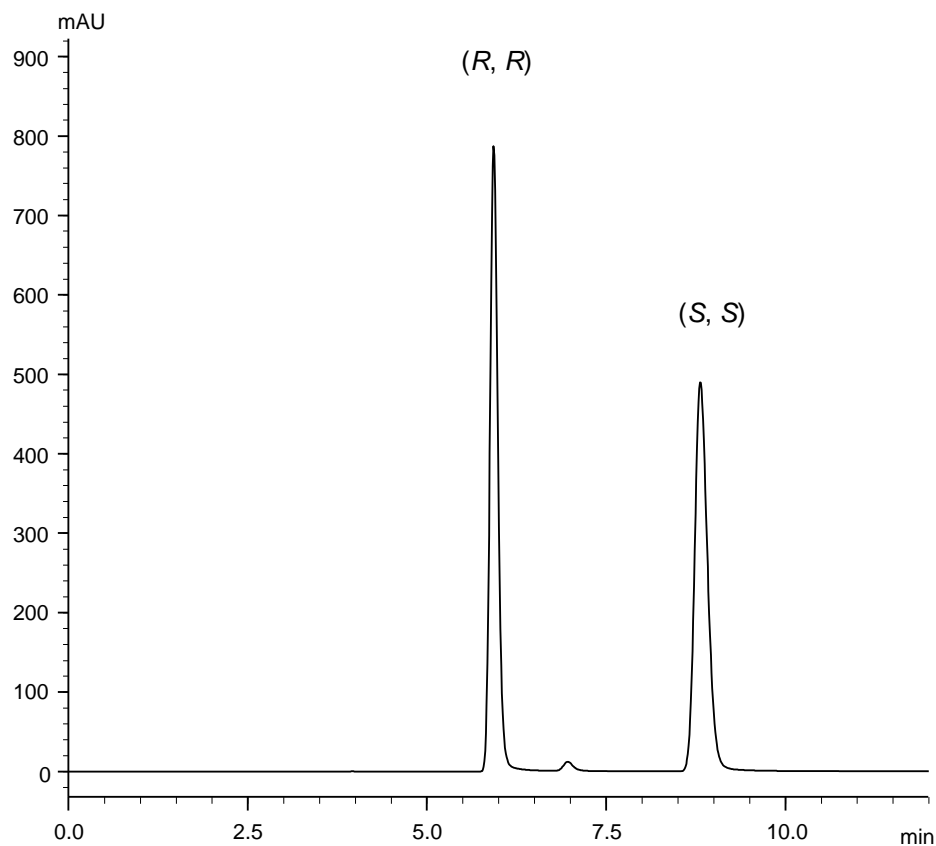

α	1.9
R_s	10.9

Column : CHIRAL ART Amylose-C Neo (5 μ m)
250 X 4.6 mm I.D.
Eluent : *n*-hexane/2-propanol (90/10)
Flow rate : 1.0 mL/min
Temperature : 25°C
Detection : UV at 285 nm
Injection : 10 μ L (0.5 mg/mL)