


抗体医薬品  
Antibody drugs

U190423A

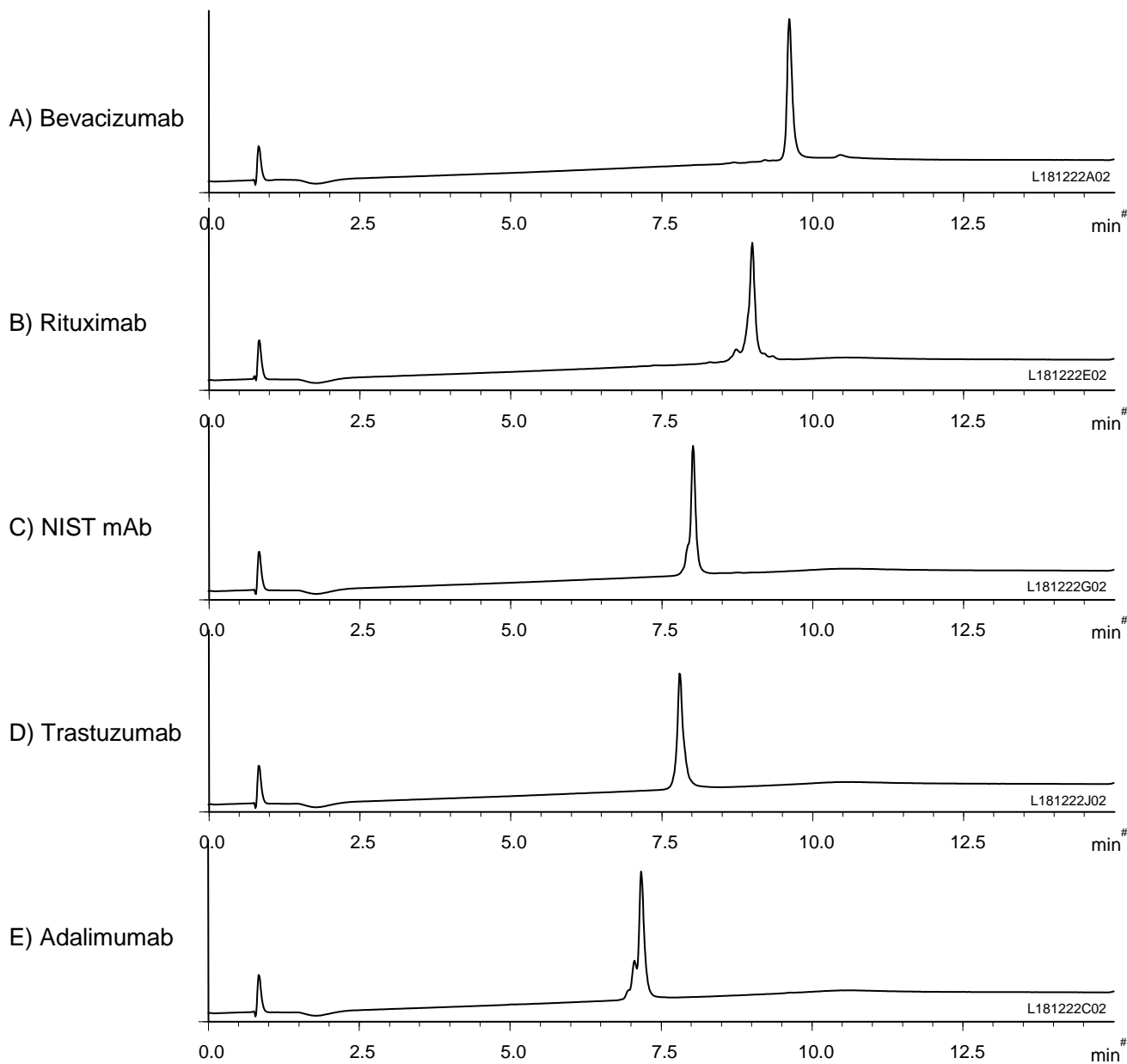

Column : BioPro HIC BF (4  $\mu$ m)  
100 X 4.6 mm I.D.

Eluent : A) 100 mM  $\text{NaH}_2\text{PO}_4$ - $\text{Na}_2\text{HPO}_4$  (pH 7.0) containing 2.0 M  $(\text{NH}_4)_2\text{SO}_4$   
B) 100 mM  $\text{NaH}_2\text{PO}_4$ - $\text{Na}_2\text{HPO}_4$  (pH 7.0)  
0-100%B (0-10 min), 100%B (10-15 min)

Flow rate : 1.0 mL/min

Temperature : 25°C

Detection : UV at 280 nm

Injection : 5  $\mu$ L (0.5 mg/mL)