



抗体医薬品
Antibody drugs

U190422B


A) Adalimumab

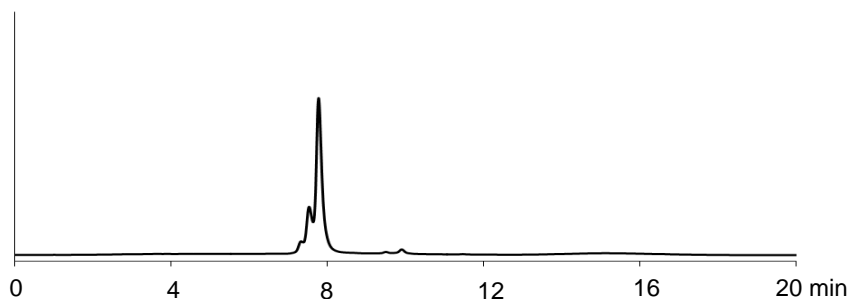


B) Brentuximab

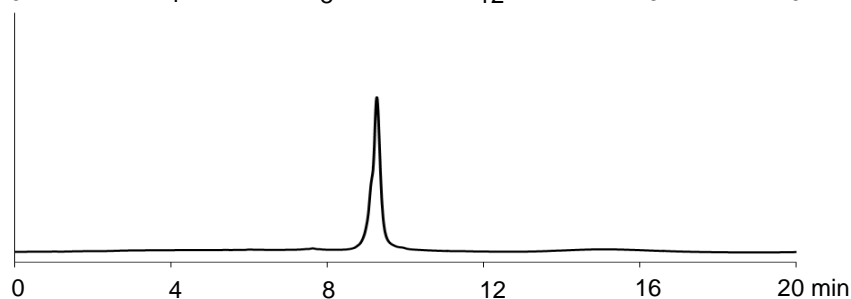

C) Pembrolizumab

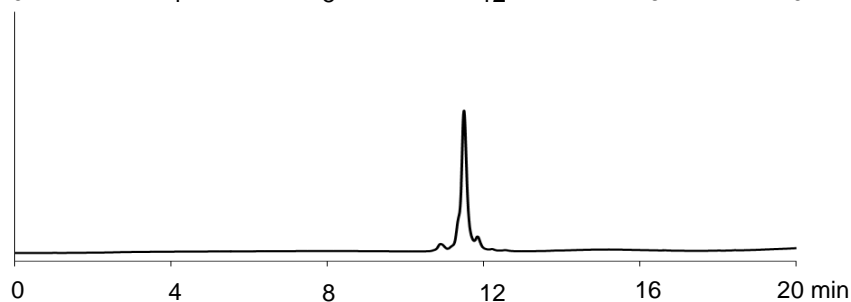

D) Reslizumab

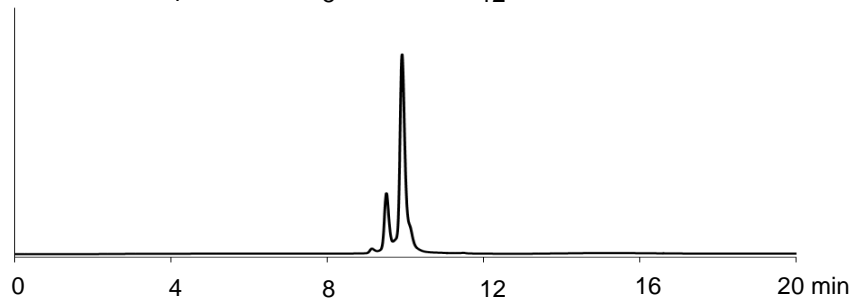

Courtesy of University of Geneva

Column	: BioPro HIC BF (4 μ m) 100 X 4.6 mm I.D.
Eluent	: A) 20 mM NaH_2PO_4 - Na_2HPO_4 (pH 7.4) containing 1.5 M $(\text{NH}_4)_2\text{SO}_4$ B) 20 mM NaH_2PO_4 - Na_2HPO_4 (pH 7.4)/2-propanol (85/15) 0-100%B (0-20 min)
Flow rate	: 1.0 mL/min
Temperature	: 20°C
Detection	: FLS at Ex 280 nm, Em 360 nm
Injection	: 3 μ L (2 mg/mL)