
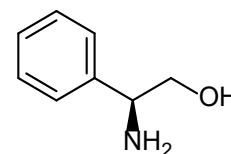



2-フェニルグリシノール

2-Phenylglycinol

E181116F

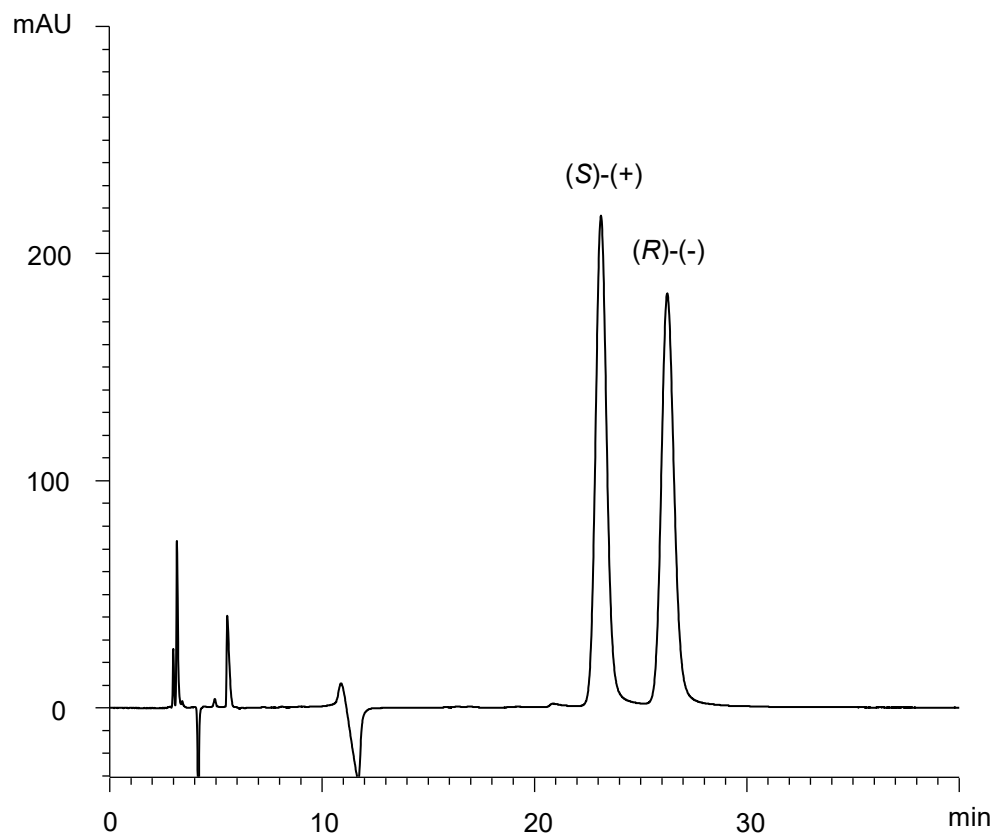

(S)-(+)

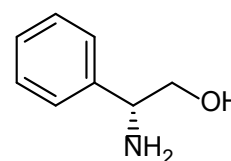

(R)-(-)

α	1.1
R_s	3.0

Column : CHIRAL ART Cellulose-SB (5 μ m)
250 X 4.6 mmI.D.

Eluent : *n*-hexane/tetrahydrofuran/trifluoroacetic acid/*n*-butylamine
(80/20/0.1/0.1)

Flow rate : 1.0 mL/min

Temperature : 25°C

Detection : UV at 215 nm

Injection : 10 μ L (0.5 mg/mL)

Sample : Supplied by Tokyo Chemical Industry Co., Ltd.