
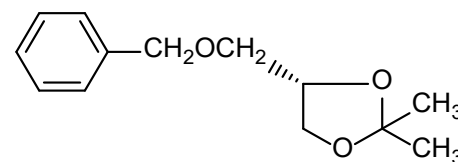



4-ベンジルオキシメチル-2,2-ジメチル-1,3-ジオキサラン

4-Benzyloxymethyl-2,2-dimethyl-1,3-dioxolane

E180517B

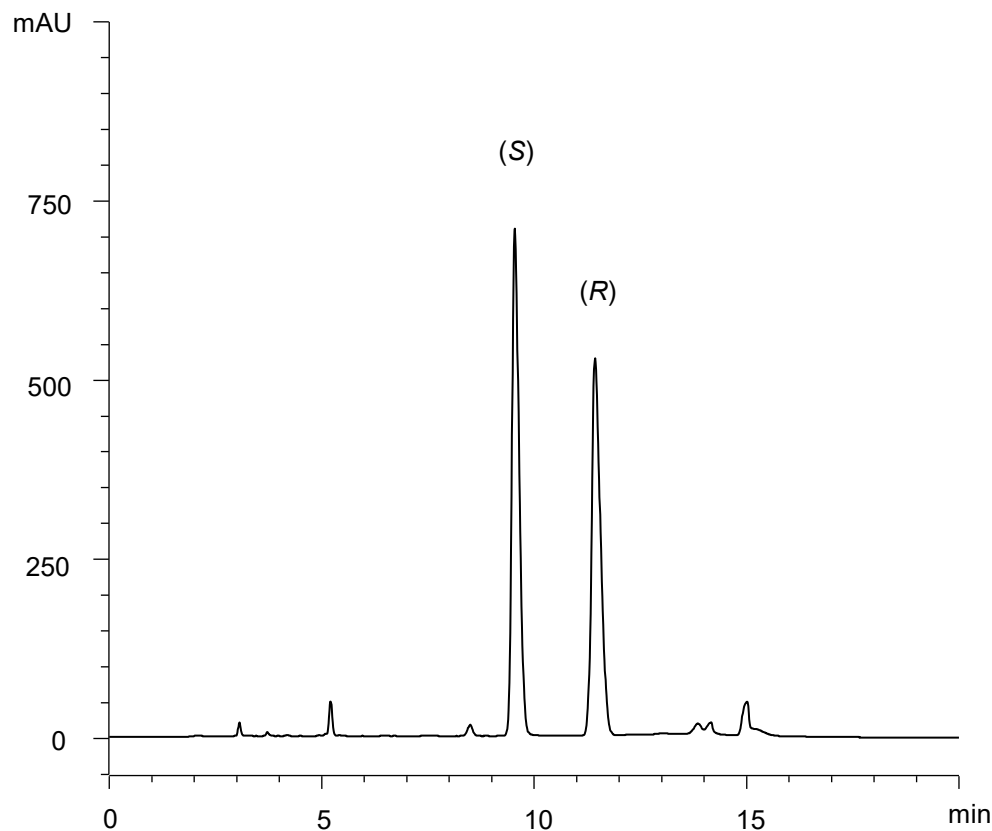

(S)

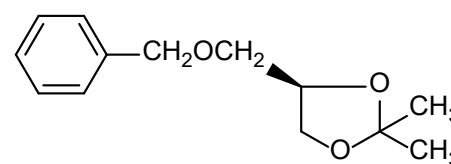

(R)

α	1.2
R_s	5.2

Column : CHIRAL ART Cellulose-SB (5 μ m)
250 X 4.6 mm I.D.

Eluent : *n*-hexane/methyl *tert*-butyl ether (90/10)

Flow rate : 1.0 mL/min

Temperature : 25°C

Detection : UV at 215 nm

Injection : 10 μ L (0.5 mg/mL)