
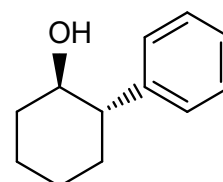



trans-2-フェニル-1-シクロヘキサノール

trans-2-Phenyl-1-cyclohexanol

E180509J

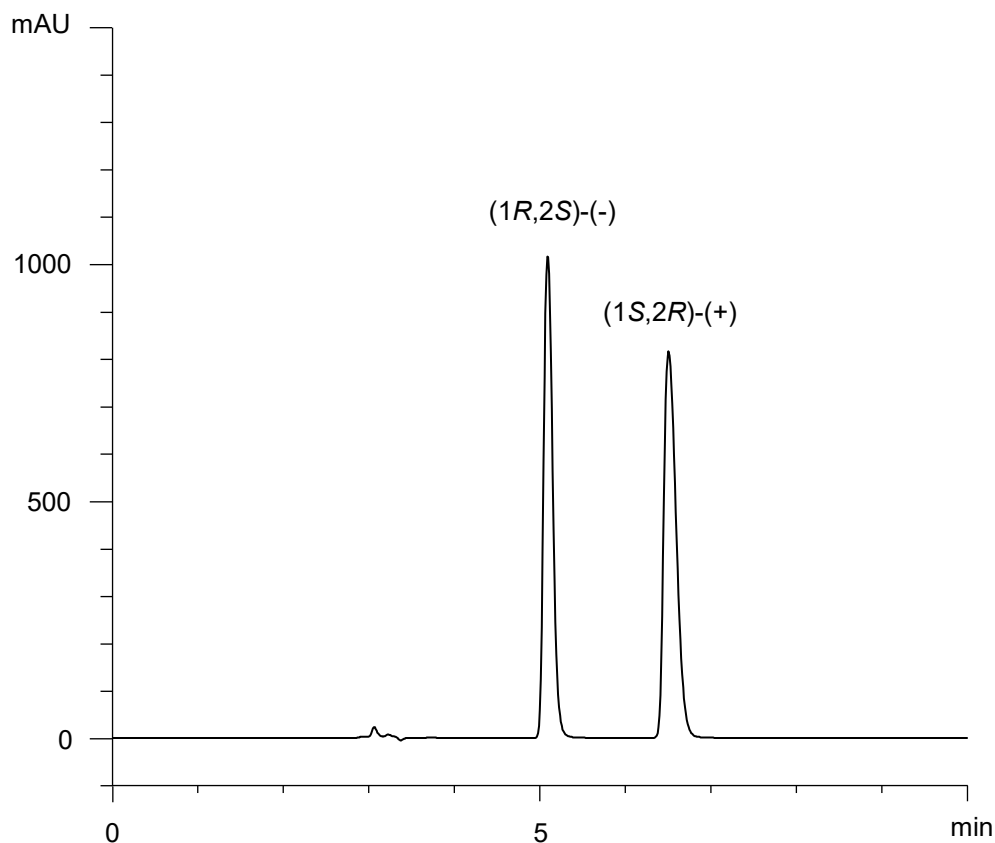

(1R,2S)-(-)

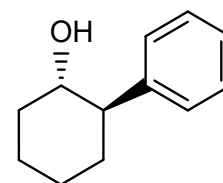

(1S,2R)-(+)

α	1.6
R_s	5.8

Column : CHIRAL ART Amylose-SA (5 μ m)
250 X 4.6 mm I.D.

Eluent : *n*-hexane/ethanol (70/30)

Flow rate : 1.0 mL/min

Temperature : 25°C

Detection : UV at 215 nm

Injection : 10 μ L (0.5 mg/mL)

Sample : Supplied by Tokyo Chemical Industry Co., Ltd.