
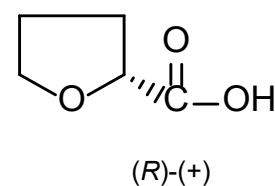
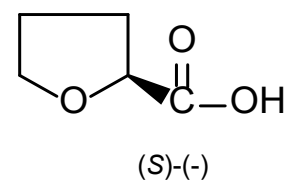


テトラヒドロフラン-2-カルボン酸
Tetrahydrofuran-2-carboxylic Acid

E171109J

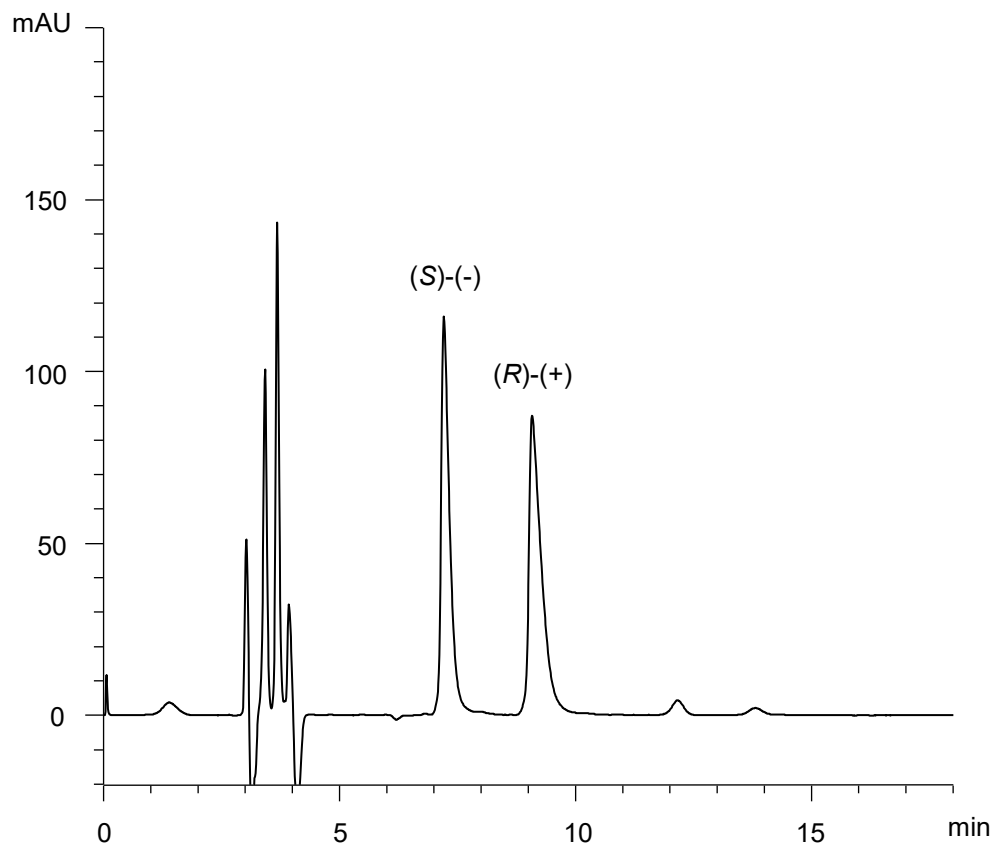

α	1.4
R_s	4.5

Column : CHIRAL ART Cellulose-SC (5 μ m)
250 X 4.6 mm I.D.
Eluent : *n*-hexane/2-propanol/trifluoroacetic acid (80/20/0.1)
Flow rate : 1.0 mL/min
Temperature : 25°C
Detection : UV at 215 nm
Injection : 10 μ L (5 mg/mL)