



モノクローナル抗体 (アダリムマブ)  
Monoclonal antibody (Adalimumab)

V180731D

Mobile phase containing 0.1% TFA

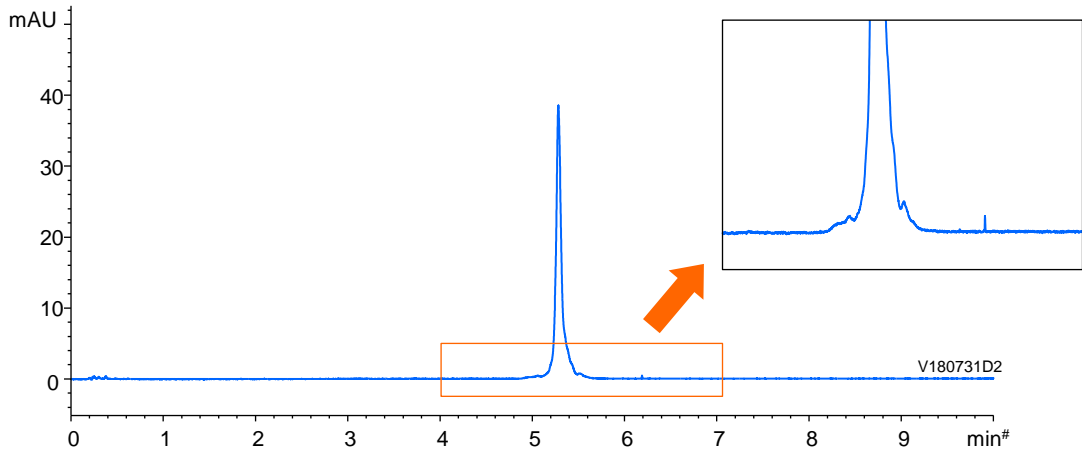

Mobile phase containing 0.1% formic acid

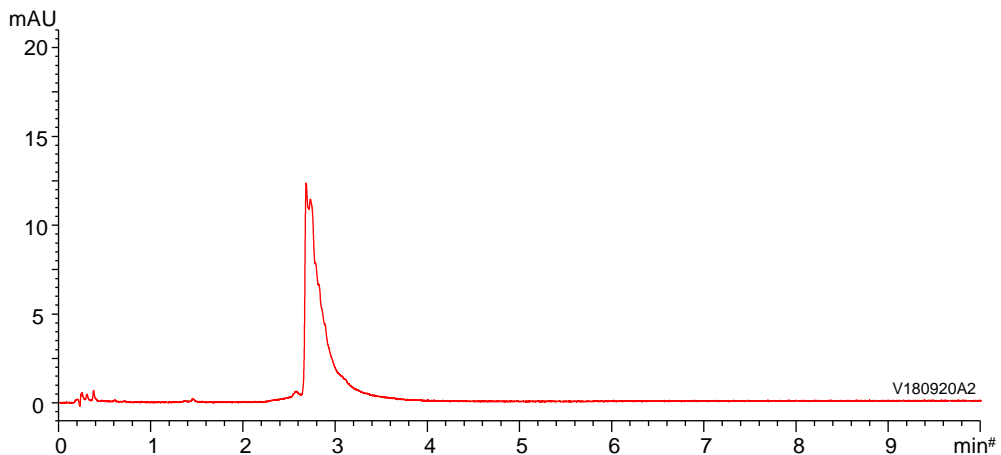

Column : YMC-Triart Bio C4 (1.9  $\mu$ m, 30 nm)  
50 X 2.1 mmI.D.

Eluent : A) water/TFA or formic acid (100/0.1)  
B) acetonitrile/TFA or formic acid (100/0.1)  
25-45%B (0-10 min), 90%B (10-12.5 min)

Flow rate : 0.4 mL/min

Temperature : 80°C

Detection : UV at 280 nm

Injection : 2  $\mu$ L (0.5 mg/mL)