



モノクローナル抗体 (ベバシズマブ)
Monoclonal antibody (Bevacizumab)

V180731E

Mobile phase containing 0.1% TFA

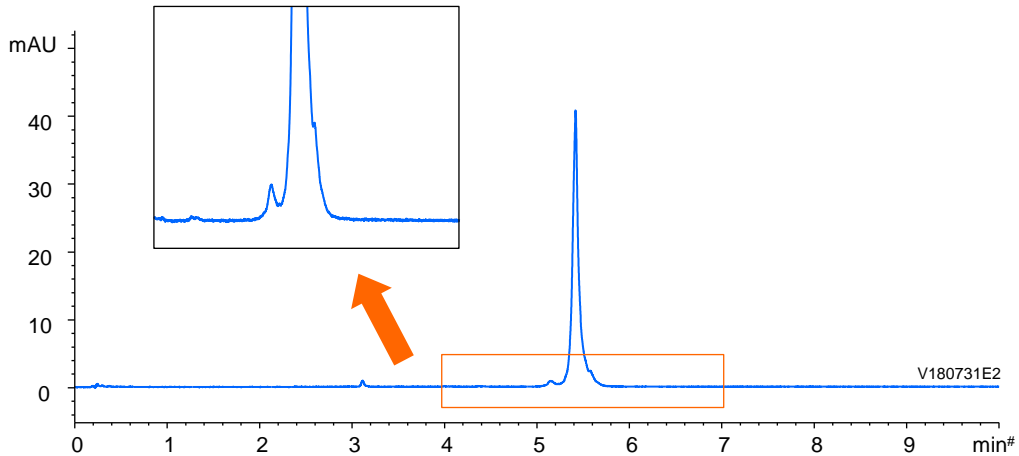

Mobile phase containing 0.1% formic acid

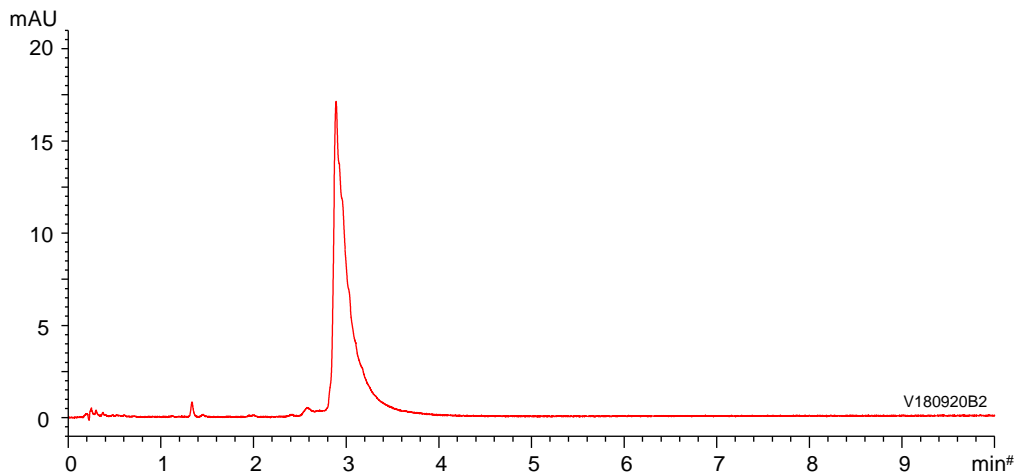

Column : YMC-Triart Bio C4 (1.9 μ m, 30 nm)
50 X 2.1 mmI.D.

Eluent : A) water/TFA or formic acid (100/0.1)
B) acetonitrile/TFA or formic acid (100/0.1)
25-45%B (0-10 min), 90%B (10-12.5 min)

Flow rate : 0.4 mL/min

Temperature : 80°C

Detection : UV at 280 nm

Injection : 2 μ L (0.5 mg/mL)