


# Application Note



Paclobutrazol

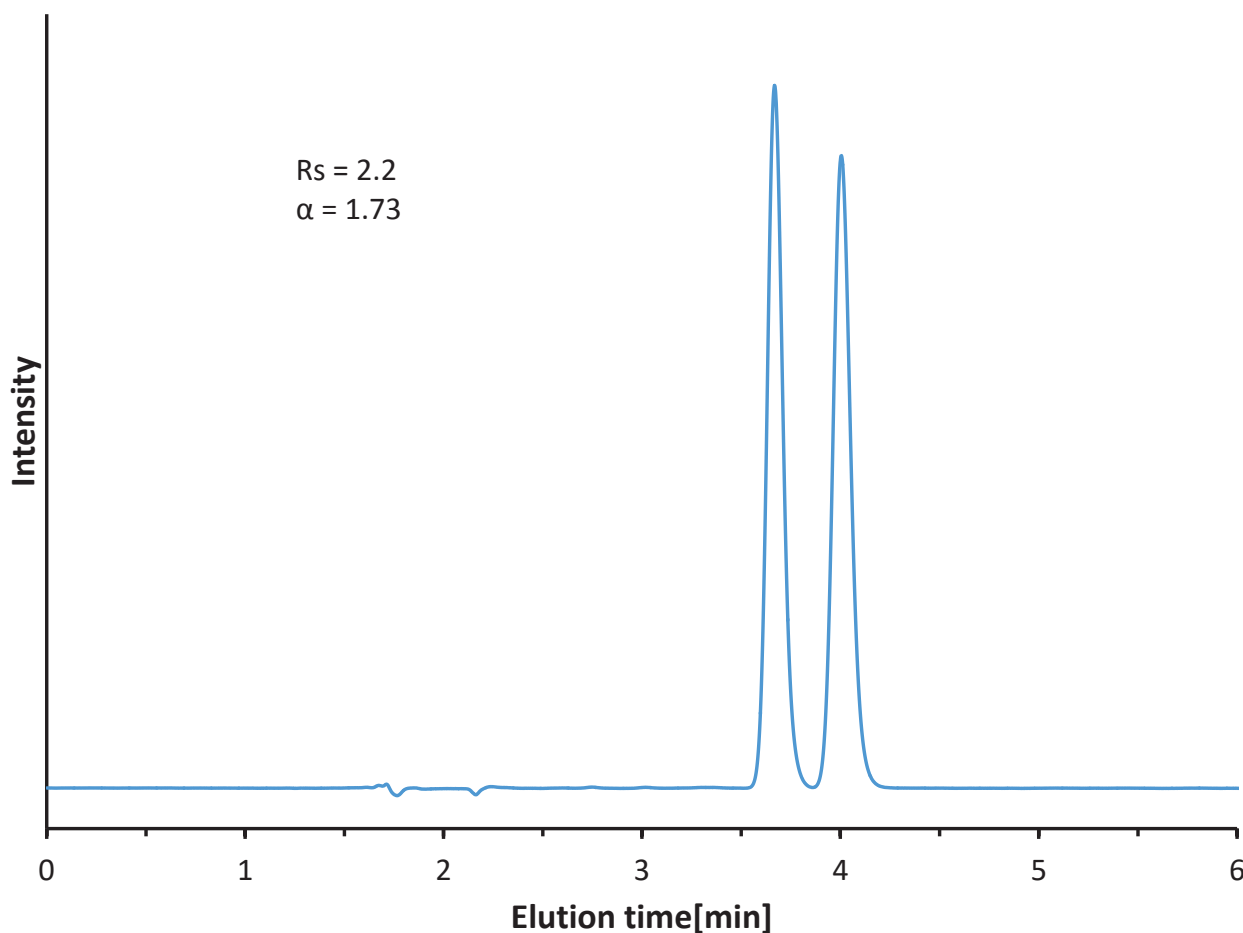

Column: CHIRAL ART Cellulose-SZ (5  $\mu$ m) 150 x 4.6 mm ID  
Part No.: KSZ99S05-1546WT  
Eluent: 20 mM  $\text{NH}_4\text{HCO}_3$  / acetonitrile / diethylamine (50/50/0.1)  
Flow rate: 1 mL/min  
Temperature: 25  $^\circ\text{C}$   
Detection: UV at 220 nm  
Injection: 5  $\mu\text{L}$   
Sample: Paclobutrazole, 1 mg/mL in eluent