



ネアラビン
Nelarabine

U211015A

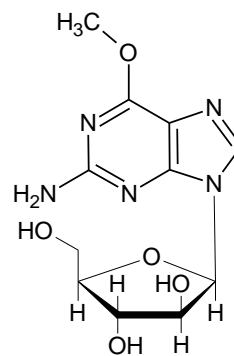

Nelarabine

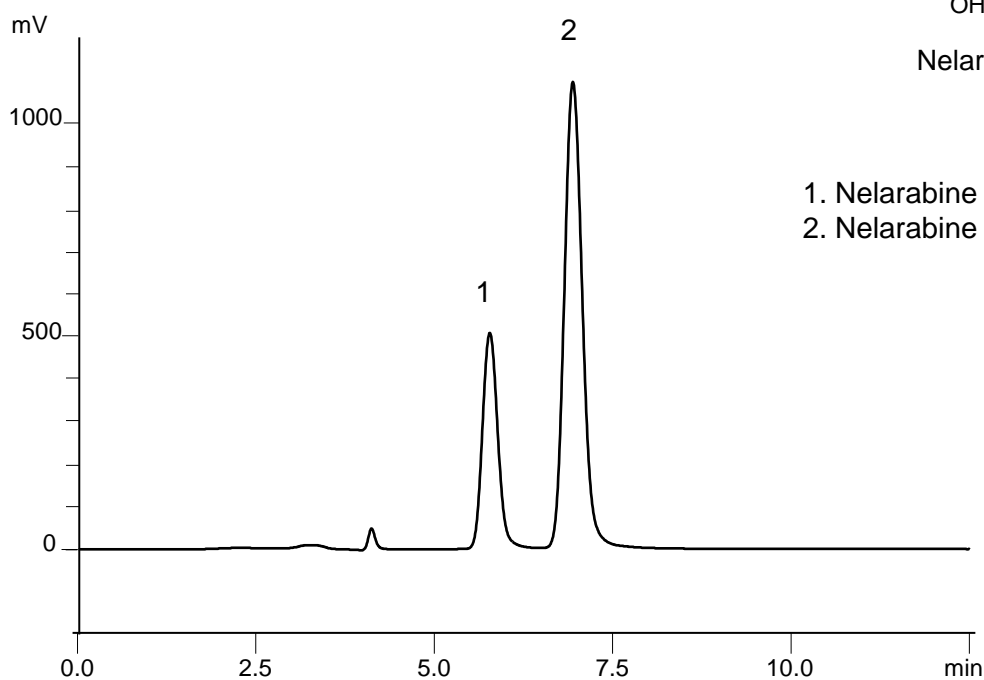

1. Nelarabine enantiomer
2. Nelarabine

Courtesy of YMC India Pvt. Ltd.

Column	: CHIRAL ART Cellulose-SC (5 μ m) 250 X 4.6 mm I.D.
Eluent	: <i>n</i> -hexane/ethanol/diethylamine (50/50/0.1)
Flow rate	: 1.0 mL/min
Temperature	: 25°C
Detection	: UV at 266 nm
Injection	: 10 μ L (1 mg/mL)