

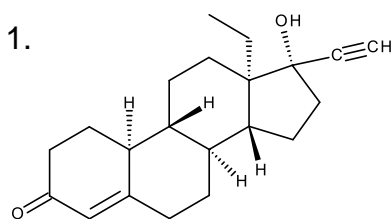


DL-ノルゲストレル  
DL-Norgestrel

AB220131A

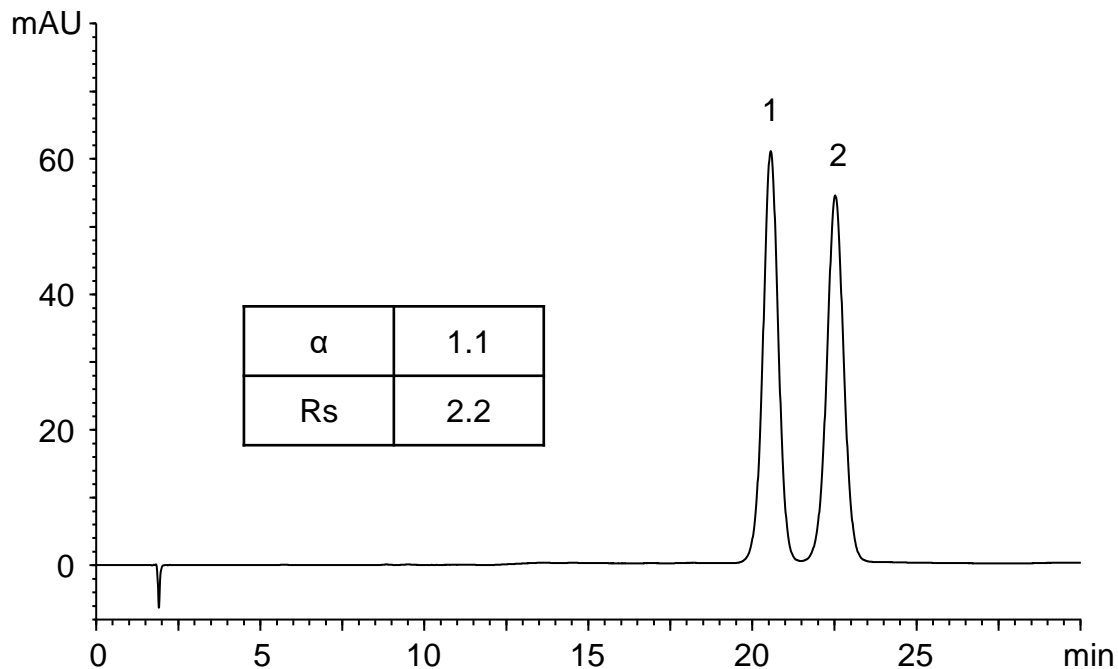

L-Norgestrel

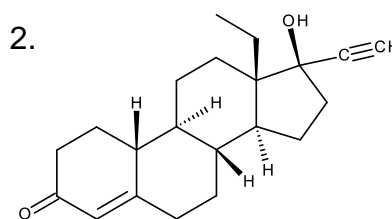

D-Norgestrel  
(Levonorgestrel)

Column : CHIRALART Cellulose-SC (3  $\mu$ m)  
150 X 2.0 mm I.D.  
Eluent : acetonitrile/water/formic acid (45/55/0.1)  
Flow rate : 0.2 mL/min  
Temperature : 25°C  
Detection : UV at 240 nm  
Injection : 1  $\mu$ L (0.5 mg/mL)