


膜透過性ペプチド  
Cell penetrating peptides

AB211202A

1. Octaarginine
2. HIV-1 Tat protein (47-57)
3. Penetratin

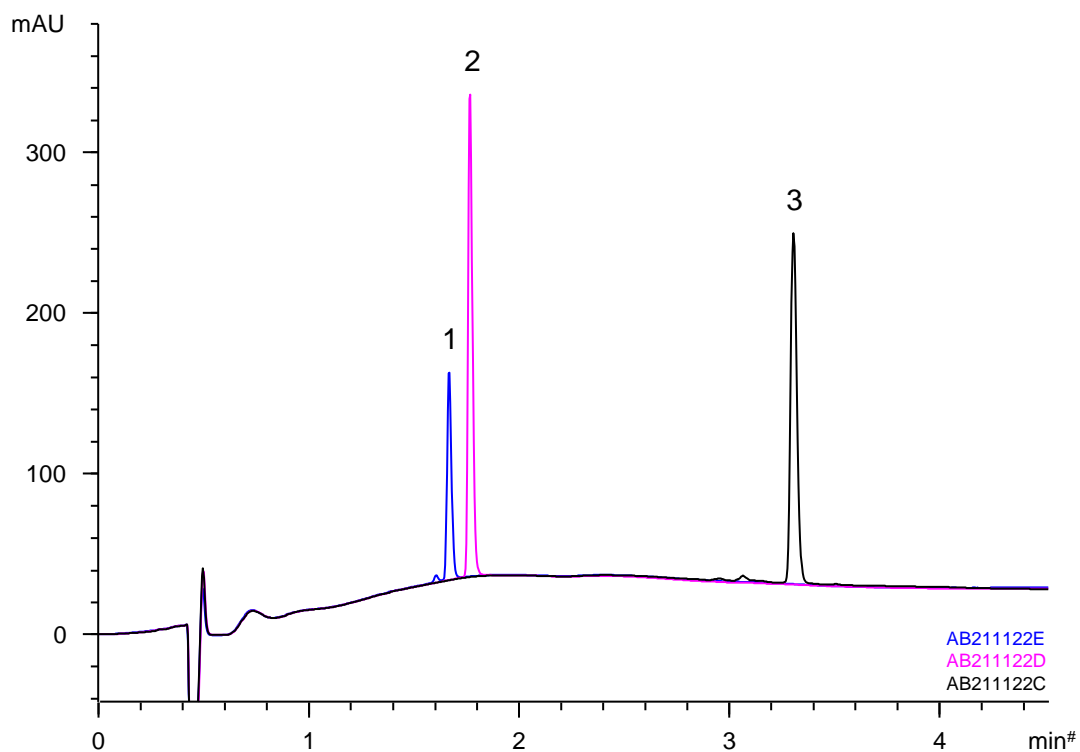

Column	: YMC-Triart C18 (1.9 $\mu$ m, 12 nm) 50 X 2.1 mm I.D.
Eluent	: A) water/TFA (100/0.1) B) acetonitrile/TFA (100/0.08) 5-25%B (0-2.5 min), 25%B(2.5-5 min)
Flow rate	: 0.4 mL/min
Temperature	: 40°C
Detection	: UV at 210 nm
Injection	: 2 $\mu$ L (0.05~0.5 mg/mL)