
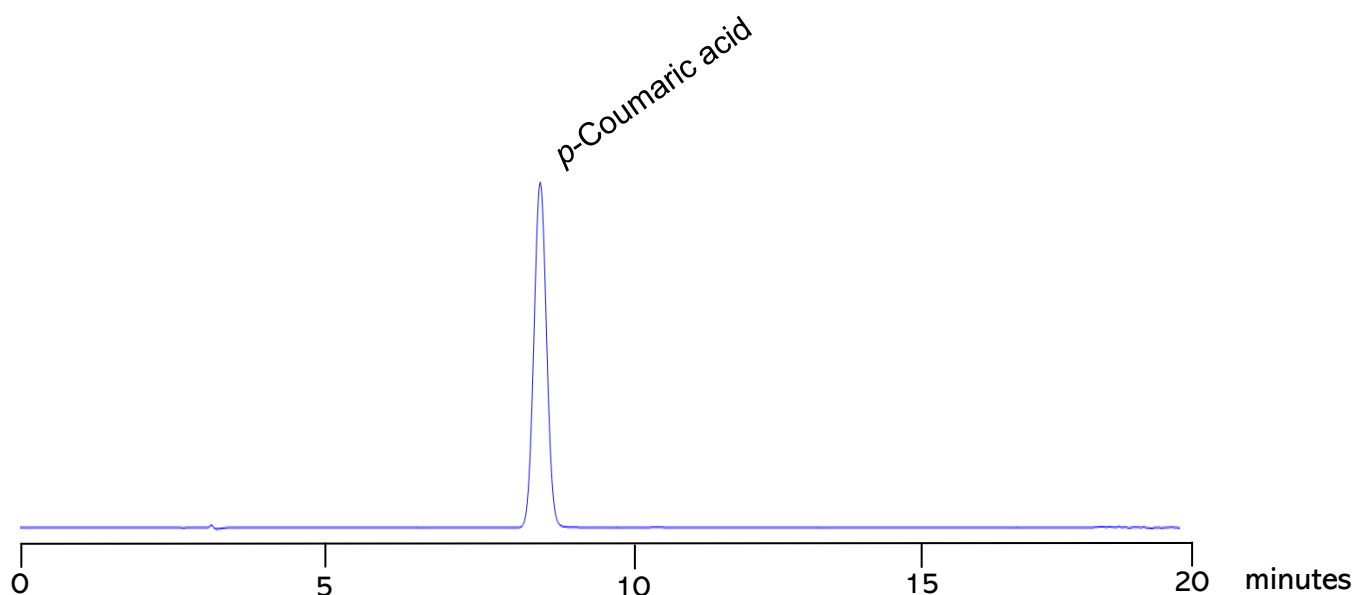


p-Coumaric acid

파라-쿠마르산

KR23051801

p-Coumaric acid

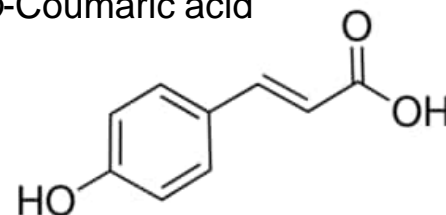

Column	: YMC Triart C18 (5 μ m, 12nm) 250 x 4.6 mm I.D,
Eluent	: A) 0.5% acetic acid B) methanol 45-50%B (0-15 min), 50-100%B (15-15.01 min), 100%B (15.01-25 min), 45%B (25.01-35 min)
Flow rate	: 1.0 mL/min
Detection	: UV at 310 nm
Temperature	: 40 $^{\circ}$ C
Injection	: 10 μ L